

Finest ENT Equipment.  
Made in Germany. Since 1950.



## BASIC PLUS Flat Top


High-end ENT treatment unit  
offering a wide range of configuration options

**otopront**<sup>®</sup>

# High-end ENT treatment unit with a wide range of configuration options

- Body available in four sizes
- Extensive configuration options
- Well-structured functional elements
- Practice-oriented workplace layout
- High-quality materials
- Reliability and durability in clinical use

## Outstanding features

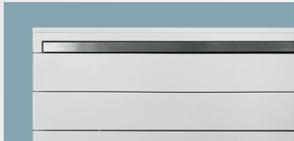

Drawer fronts are also available in high-quality mineral material with push-to-open function

# Instrument management and storage



## Instrument Trays

- Ergonomic, hygienic, and organized instrument handling
- Dust- and splash-proof storage
- Removable trays with antibacterial coating
- Adjustable stainless-steel dividers for flexible setup
- All components are mechanically reprocessable



## Drawers

- Flexible layout design:
  - Clean/unclean or storage/disposal configuration
  - Customizable to room and workflow requirements
- Push-to-open / Softclose system available



## Plateau

- Additional storage
- Easily accessible while standing
- Surrounding edge prevents items from falling off

# Suction, LED light, air pressure

## Suction Unit

- Basic pump: 58 Ltr./min., 91% vacuum  
*Alternatively:*
- Diaphragm pump: 72 Ltr./min., 91% vacuum
- Double-piston pump: 55 Ltr./min., 93% vacuum, designed for continuous operation
- All pumps are vibration-free and extremely quiet running
- Vacuum gauge
- Double overflow protection
- Automatic activation, optionally by foot switch
- Plug-in suction glass
- Use of disposable suction bags and suction hoses, no modifications required

## Suction Hose Cleaning

- Integrated rinsing of hose and cannulas with disinfectant solution

## LED Light Sources

- 1 or 2 LED light sources (CRI >95%, >1,800 lm, adjustable brightness, warm white)
- Otopront Superbright-LED (>2,600 lm)

## Spray Unit

- 0–4 bar, adjustable, with gauge
- Fine dosing
- Spray handpiece, autoclavable
- Two spray bottles with replaceable cannulas
- Otopront disposable cannulas

Customized  
solutions for  
all your needs



# BASIC PLUS Flat Top Configuration overview

## 1 Choose from 4 Dimensions



Small body  
Standard height



Small body  
Extended height



Wide body  
Standard height



Wide body  
Extended height

▼ The following configurations are available for **both widths** and **both heights** ▼

## 2 Base Functions – Suction, LED Light, Air Pressure



Suction only



Suction & LED light



Suction & air pressure  
(optional LED light)



Functions laterally reversed

## 3 Microscope & Support Options > Units also available without plateau



Small plateau & microscope



Big plateau & microscope



Swivel support & microscope



Plateau & microscope laterally reversed

## 4 Camera Mounting Options



With PES Eco HD not visibly integrated  
(optional wall mounted screen)



Small plateau, microscope with  
PES Eco HD & screen



PES PILOT HDpro on swivel support  
& microscope



PES PILOT HDpro  
on big plateau & microscope


## BASIC PLUS Flat Top Additional features



High-quality mineral material surface with protecting edge flush with the back wall



Available in many patterns, surface designs, and in almost all RAL colors

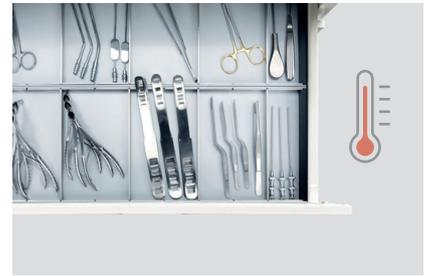

Heatable top drawer for body-warm instruments



Lighting for work surfaces, mounted under the plateau




Lateral headlight suspension



Metal holder for dispensers of cotton wool, tongue holder tissues and similar items, mounted on the plateau, freely movable



Extra silent & extra powerful surgery pump with boost function

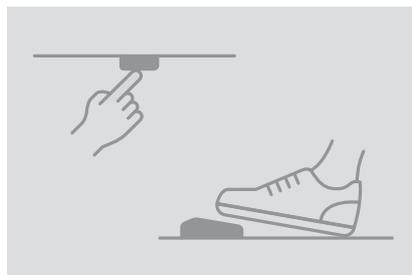

Suction on/off via button and foot switch



Integrated otoscope holder



Drawers are available in various heights and in different quantities



Mobile thanks to 360° swivel casters



Disposable spray cannulas for compressed air

# Enhance your ENT workstation with these compatible solutions

## SIT MAX one

ENT examination chair



- Spacious seat with durable, high-quality steel frame
- Height adjustment up to 400 mm (15.7 inches)  
Available via push button or foot pedal
- Mobility with electric lifting casters
- Electrically or manually adjustable backrest
- Programmable positions with an automatic reset-to-zero function
- Surgical headrest for optimal head positioning, freely adjustable (optional)
- Available in five upholstery colors: Black, Tin, Ocean Blue, Sand, Cognac
- Upholstery stitching can be customized with clinic logos, symbols, or other desired designs


## PES PILOT HDpro

Full-HD video endoscopy system




- Seamless integration directly at the Otopront ENT workstation
- Full-HD Camera (optional 4K)
- LED or stroboscopy light (optional)
- Easy-to-use, compact design
- Complete data storage & patient management

# Dimensions **without** plateau

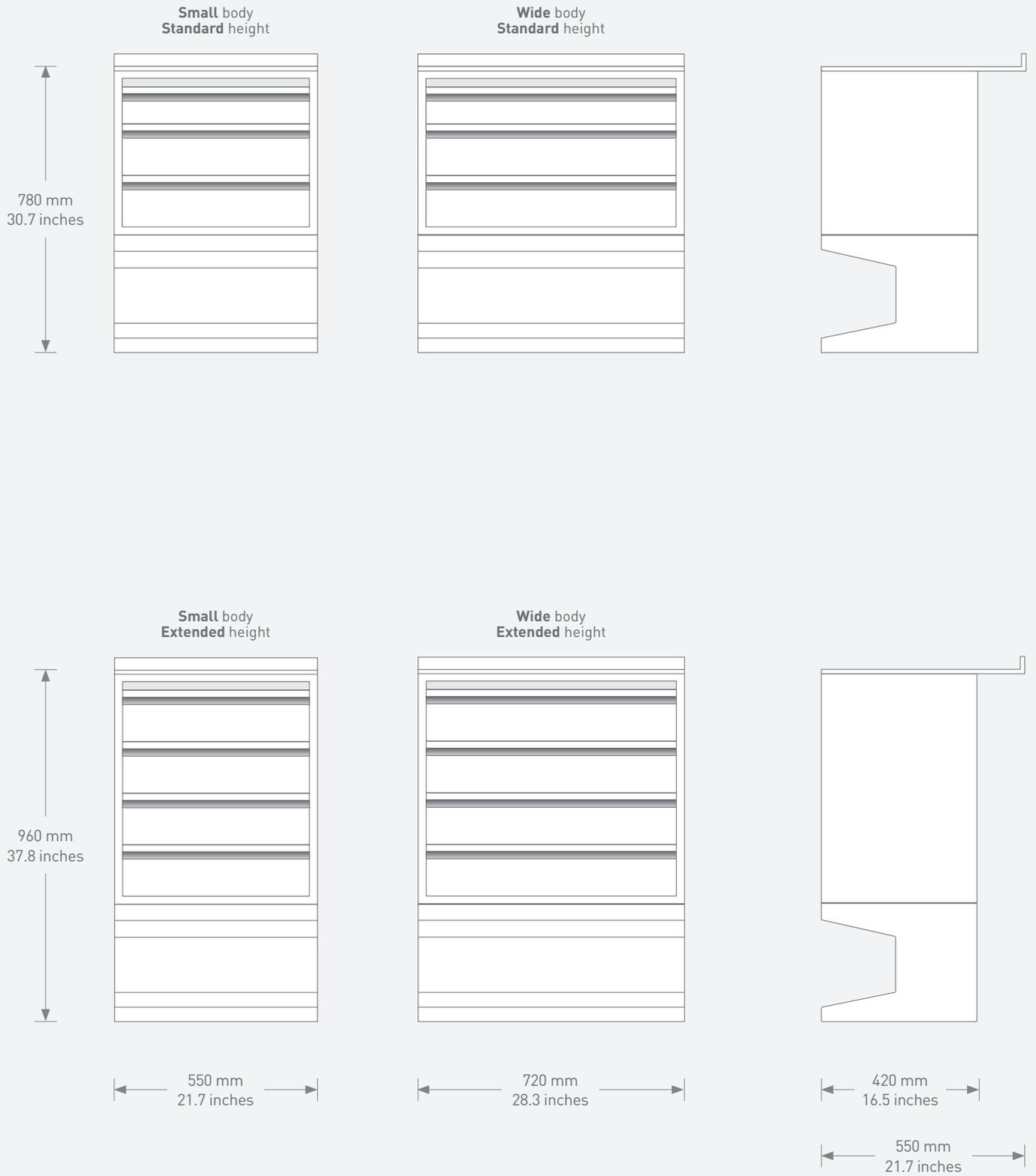

# Dimensions **with** plateau

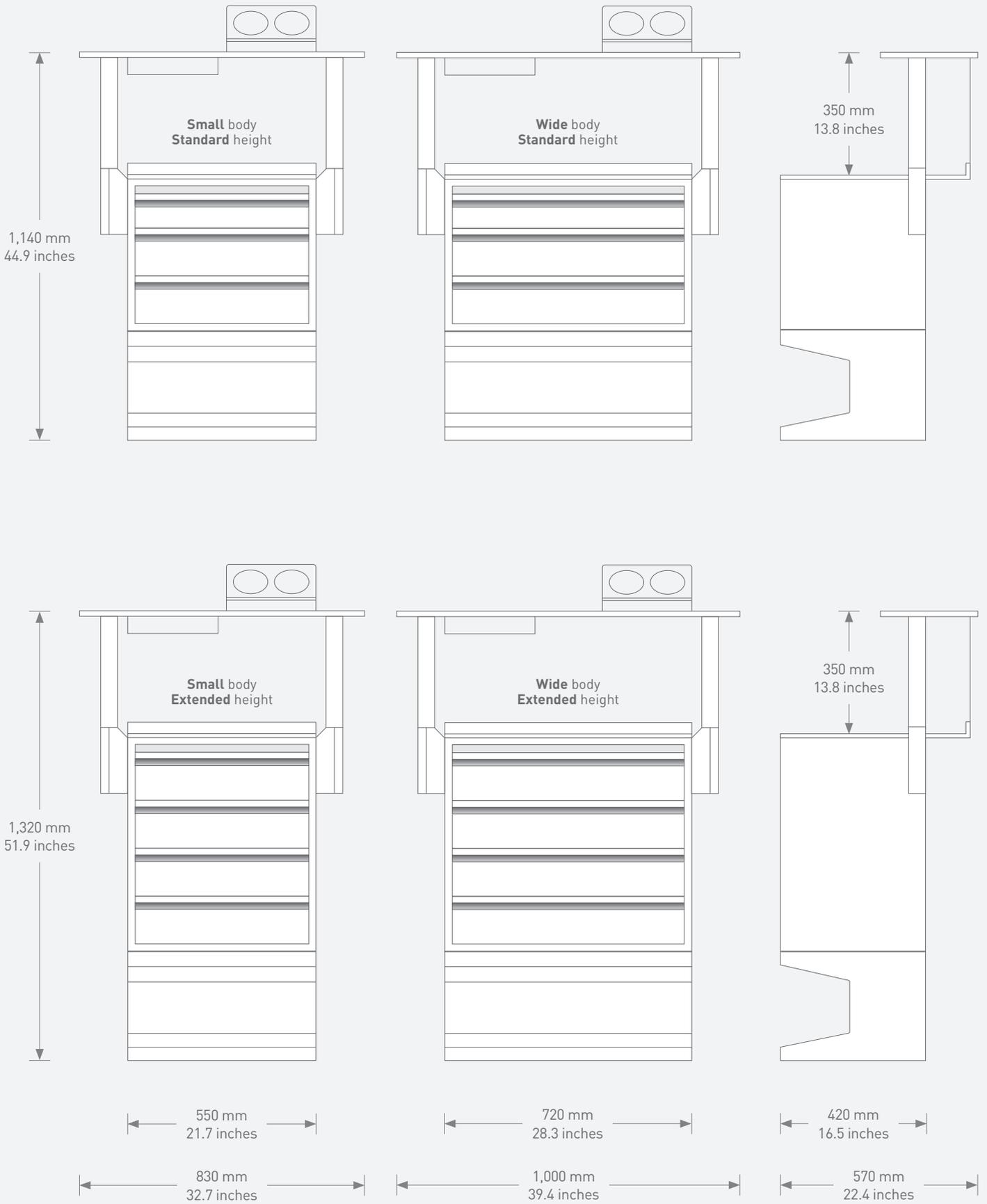

# Technical data

## Body

- Available in 4 dimensions
- High-quality steel housing
- White mineral surface plate
- Instrument drawers (with aluminum trays)
- Big storage drawers
- Drawer fronts in high-quality mineral material Corian® (optional)

## Base Functions

### Powerful Suction

- Silent, vibration-free operation
- 91% max. vacuum, 58 Ltr./min.
- Easy-to-change disposable bags
- Side-mounted disposable suction canister

### LED Light

- Otopront super bright LED light source
- CRI >95%, >1,800 lumens
- Adjustable brightness, warm white

### Air Pressure

- Compressed air system
- Fine dosing spray handpiece, flow limiter and pressure gauge
- 0–4 bar, adjustable, with gauge
- Two spray bottles with replaceable cannulas

## Microscope

- *Diagnostic Microscope ZEISS OPMI pico®*
- Mounted at the ENT Unit (incl. rod)
- 5-step magnification changer
- Straight tube f= 170 mm
- Two 10x widefield eyepieces
- Objective lens f= 250 mm
- T-handle at front

## Support Options

- Small plateau
- Big plateau
- Swivel support

BASIC PLUS Flat Top is only available for the US market.  
All designs and specifications subject to change without notice.  
All pictures shown are for illustration purpose only.  
Actual product may vary due to product enhancement.



350 mm  
13.8 inches

960 mm  
37.8 inches

420 mm  
16.5 inches

720 mm  
28.3 inches

550 mm  
21.7 inches\*

\* Body is also available in a width of 550 mm / 21.7 inches



Discover more solutions on [www.otopront.com](http://www.otopront.com) and [www.otopront.us](http://www.otopront.us)

**Happersberger otopront GmbH**

Medical Equipment Manufacturer  
65326 Hohenstein • P.O. Box • Germany  
Tel. +49 6120-9217-0 • Fax +49 6120-9217-60  
[www.otopront.us](http://www.otopront.us) • [info@otopront.com](mailto:info@otopront.com)

45968US

**otopront®**